




DERMASUCTION™

Pore Cleaning Device

DELUXE RECHARGEABLE EDITION

IMPORTANT SAFETY INSTRUCTIONS. READ AND FOLLOW ALL WARNINGS AND INSTRUCTIONS BEFORE USING THIS PRODUCT. INJURY CAN RESULT FROM IMPROPER USE. SAVE THESE INSTRUCTIONS.

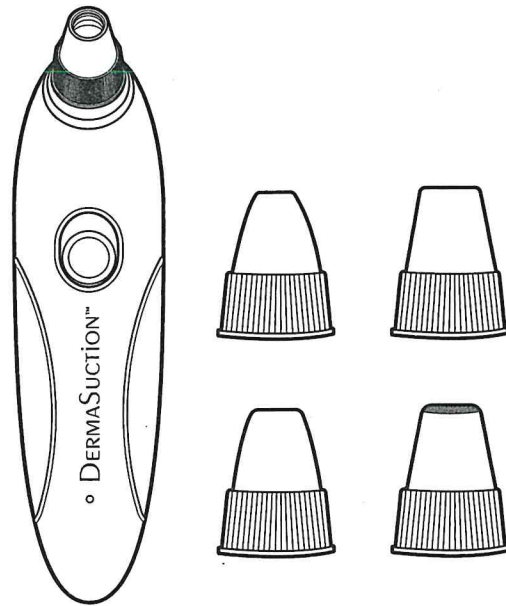

PRODUCT MAY VARY FROM IMAGES SHOWN

TABLE OF CONTENTS

Warnings.....2
Product Components.....3
How To Charge.....4
How To Change The Suction Heads.....4
How To Use.....5
How To Clean.....6
FCC Compliance.....6



WARNING

- Before use, talk to your doctor if you have any medical issues or concerns.
- This product is not for use by children.
- Do not wear necklaces, earrings or any facial jewelry while using this product.
- Necklaces, earrings and other types of jewelry should be removed before using this product.
- This product is not for use on areas of the body that have undergone any sort of surgery or have any scars.
- If you experience pain while using this product, discontinue use immediately.
- This product is not water proof and should not be used in the shower or any situation that involves water.
- Do not store this product in the bathroom or a room that has high moisture content.
- Do not operate in the presence of explosives and/or other flammable fumes or liquids.
- Do not store or use this product near sources of heat or fire.




WARNING



STORE OUT OF REACH OF CHILDREN.
THIS PRODUCT IS NOT A TOY.

PRODUCT MAY VARY FROM IMAGES SHOWN

PRODUCT COMPONENTS

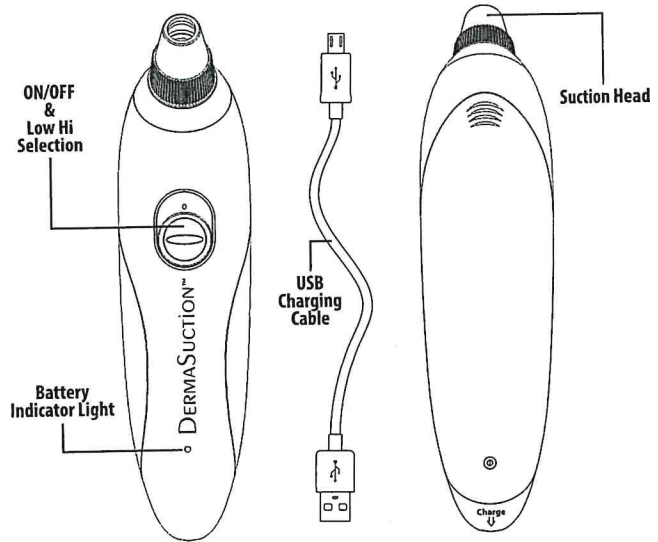

Small Head



Gentle suction, perfect for everyday use. Good for use on sensitive skin and nose area.

Large Head



Strong suction, perfect for blackheads.
Note: Not for use around the eye area.

Oval Head



Oval opening for getting into corners of nose and mouth, and for sensitive skin around eyes.

Sonic Head




Helps exfoliate dead skin while pulling dirt & oil from skin.
Note:
For Dry skin- once a week
For Oily/Mixed skin- twice a week
No More than 5 minutes a day

HOW TO CHARGE


1. Turn the Derasuction™ on to check the battery. Battery indicator light is solid green when fully charged. When indicator light is solid red, the battery is low.
2. To charge, plug the micro-USB connector end of the supplied USB cable into the charging port of the Derasuction™. (Figure A)
3. While charging the red light will be solid, indicating the Derasuction™ is charging properly. When charging is complete the light will revert back to solid green. (Figure B)

Figure A

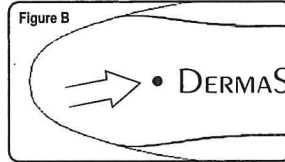


NOTE: MAKE SURE TO REMOVE THE USB CABLE FROM BOTH THE CHARGING PORT OF THE DERASUCTION™ AND THE USB PORT OF THE POWER SOURCE AFTER THE DERASUCTION™ IS FULLY CHARGED. OVER-CHARGING THE BATTERY CAN CAUSE POOR BATTERY LIFE.

HOW TO CHANGE THE SUCTION HEADS

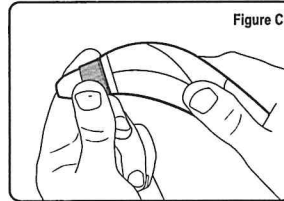


1. Place one hand on the suction head and your other hand on the unit. (Figure C)

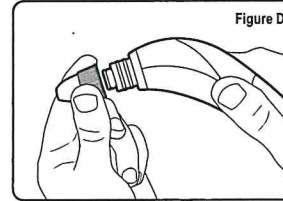


2. Firmly pull off the suction head and reinstall your desired suction head in the same manner. (Figure D)


HOW TO USE

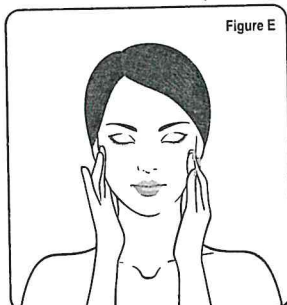

1. Start by opening up your pores by washing your face with warm water, steam or a hot wet towel. (Figure E) Make sure to dry face and hands completely before moving onto the next step.

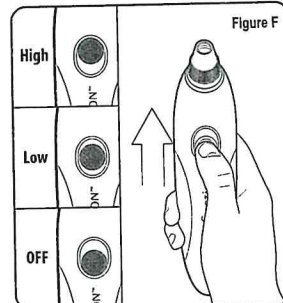

2. Select your desired suction head (see page 4 for how to change the suction head). DermaSuction™ is turned on by sliding the ON/OFF switch up. The middle position is low suction and the top position is high suction. (Figure F)

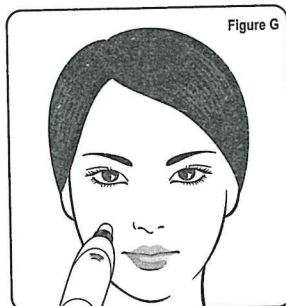


3. Apply the DermaSuction™ to your face. (Figure G) Slowly move the DermaSuction™ around your face, not keeping it in one place for more than 2 seconds. For best results, focus on your problem areas. *Note: DO NOT move back and forth in same direction.*

Helpful Tips

- After using DermaSuction™, make sure to apply moisturizer to the treated area.
- Before and after use, do not use any facial scrubs or other exfoliating products.
- Do not apply any downward force during use, the suction capabilities of DermaSuction™ are enough to clean your pores efficiently.
- Slight redness is normal and will disappear after approximately 30 minutes.

HOW TO CLEAN

1. After each use, clean the suction head by running under warm water. (Figure H) Allow the suction head to air dry. Do not attempt to clean the main unit by running under water. To clean the main unit, wipe down with a soft damp cloth.

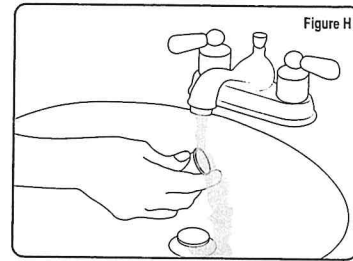

FCC COMPLIANCE

This device complies with Part 15 of the FCC Rules.

Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the Instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver connected.
- Consult the dealer or an experienced radio/TV technician for help.

Manufactured for
TELEBrands.
Fairfield, NJ 07004
©2018 Telebrands
www.telebrands.com

Made in China

LIMITATION OF LIABILITY

Liability limited to the purchase price of this product. Telebrands Corp. shall not be liable for any incidental or consequential damages or any implied warranty on this product. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.